



# NITINOL STONECATCHER



NiTi 4xDØ - KØ - Form - x

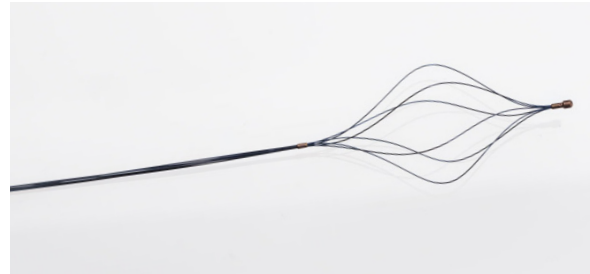

NiTi 4 (6) - DØ - KØ - Form

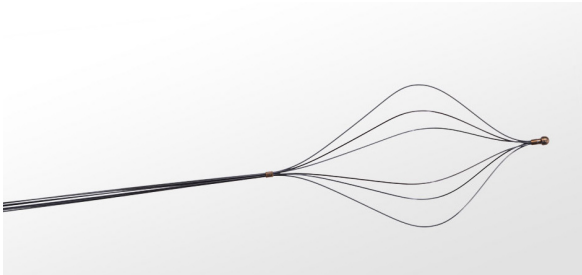

NiTi 4 (6) - DØ - KØ - Form

## Explanations

DØ	diameter of wire
KØ	diameter of basket (15, 20, 25, 30, 35 mm)
Form	GF: twisted form DF: diamond form
x	Litho: lithotripsy capable Fü: guide wire capable (only KØ 20, 25, 30, 35)

## Construction

**Basket** 4 or 6 Nitinol wires combined into a stainless steel tip.

**Pull rope** made of a stainless steel strand or continuous Nitinol wire (for lithotripsy capable basket) and stainless steel pull rod.

## Customer requirements

In addition to our standards we also manufacture stonecatchers according to your specifications and wishes.

