


# TIPLESS NITINOL STONE CATCHER

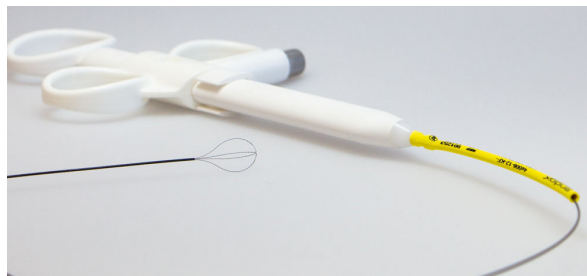

## Construction

The tipless stone catcher has 4 Nitinol wires. The nitinol pullwire is covered by a Polyimide sheath and connected to a 3-finger handle which can open, close and rotate the basket. The stone catcher is available in sizes 1.9 to 3.0 French.

- no additional components on the tip
- high tensile strength
- high expansion force
- good torque properties
- proportional opening

## Explanations

Length	90 / 120 cm
CupØ	6 / 10 / 12 / 14 mm
Material	Nitinol
Construction	4 wires
SchaftØ	1.9 / 2.6 / 3.0 Fr
Material	Polyimide / PEEK
Handle	3-finger

## Customer requirements

In addition to our standards we also manufacture tipless nitinol stone catchers according to your specifications and wishes.

